



GIMA

BORSA GHIACCIO
ICE BAG
VESSIE À GLACE
EISTASCHE
BOLSA HIELO
BORSA DE GÉLO
ΘΕΡΜΟΦΟΡΑ ΠΑΓΟΥ

MANUALE D'USO E MANUTENZIONE
USE AND MAINTENANCE BOOK
INSTRUCTIONS DE FONCTIONNEMENT ET ENTRETIEN
BETRIEBS UND WARTUNGS ANWEISUNGEN
MANUAL DE USO Y MANTENIMIENTO
MANUAL DE USO E MANUTENÇÃO
ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΗΣ ΚΑΙ ΣΥΝΤΗΡΗΣΗΣ

ATTENZIONE: Gli operatori devono leggere e capire completamente questo manuale prima di utilizzare il prodotto.

ATTENTION: The operators must carefully read and completely understand the present manual before using the product.

AVIS: Les opérateurs doivent lire et bien comprendre ce manuel avant d'utiliser le produit.

ACHTUNG: Die Bediener müssen vorher dieses Handbuch gelesen und verstanden haben, bevor sie das Produkt benutzen.

ATENCIÓN: Los operadores tienen que leer y entender completamente este manual antes de utilizar el producto.

ATENÇÃO: Os operadores devem ler e entender completamente este manual antes de usar o produto.

ΠΡΟΣΟΧΗ: Οι χειριστές αυτού του προϊόντος πρέπει να διαβάσουν και να καταλάβουν πλήρως τις οδηγίες του εγχειριδίου πριν από την χρήση του.




GIMA Spa - Via Marconi, 1 - 20060 Gessate (MI) - Italia

ITALIA: Tel. 199 400 401 (8 linee r.a.) - Fax 199 400 403

E-mail: gima@gimaitaly.com - www.gimaitaly.com

INTERNATIONAL: Tel. ++39 02 953854209 - Fax ++39 02 95380056

E-mail: export@gimaitaly.com - www.gimaitaly.com

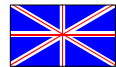

INTENDED USE

Before using it, ALWAYS ask for your doctor's opinion.

BEFORE EACH USE and from time to time apply to your doctor, who will show you possible collateral effects due to the specific use of the Device.

The present device aims at easing symptoms and/or dysfunctions generated on the human body by diseases; The ice to be used inside the bag can be produced using only drinkable water.

It can only be applied to external healthy parts of the body.

Other uses are considered improper and therefore FORBIDDEN:

- Do not use it in contact with injured parts or on parts inside the human body.
- Do not use it as invasive tool or together with invasive tools.
- Do not use together with or as part of other medical instruments, or together with drugs or substances considered improper.
- Do not modify or alter its components and/or features.
- Do not blow into the device and do not inflate it.

PRECAUTIONS AND WARNINGS

- STRICTLY OBSERVE THE INTENDED USE OF THE DEVICE.
- THE DEVICE CANNOT BE USED BY CHILDREN, DISABLE OR IMPAIRED PEOPLE WITHOUT SUPERVISION.
- IN CASE OF CONTEMPORARY USE OF SIMILAR DEVICES DO NOT INVERT DETACHABLE PARTS AND/OR ACCESSORIES.
- AFTER THE BAG IS FILLED, MAKE THE EXCEEDING AIR LEAK BEFORE PLUGGING IT. CLOSE THE BAG AND CHECK THE LOADING INLET TIGHTNESS.
- DO NOT USE OR LAY THE DEVICE ON/WITH: HOUSEHOLD ELECTRIC APPLIANCES, HOT OR POINTED SURFACES.
- BEFORE USING THE DEVICE ALWAYS MAKE SURE THAT ALL COMPONENTS ARE INTEGRAL AND NOT DAMAGED (FOR EXAMPLE CRACKING, DENTING, TEARING), THEIR COLOUR IS UNCHANGED AND THEY ARE FREE FROM SURFACE MOULDS, ETC.

FEATURES

Device class I (93/42/CEE Directive)	Device service life (normal use) See label	Capacity -
Disposable device NO	Custom-sized device NO	Device delivered STERILE NO

For a longer-lasting service life, it is advisable to store it in the original packing and in dry environment with temperatures ranging from +1 °C to + 40 °C, far from heat sources (sun or adverse weather conditions), sheltered from dust and dangerous materials and far from the children's reach.

MAINTENANCE AND DISPOSAL

Before using the device for the FIRST TIME and each time it is necessary, wash the interior of the bag with warm drinkable water and the external surface with warm water and neutral soap; do not use alcohol, solvents or acids. Do not use brushes, but only soft sponges so that not to scratch the rubber surface. After the use, completely empty the bag and make it drain/dry.

In case of replacement or breakages, dispose it according to the local legislation on the subject of waste disposal. No toxic substances are present.



GUARANTEE

Thank you for having acquired our product.

This product meets the most stringent requirements regarding the selection of quality manufacturing materials and also the final control. The product has a 12-month warranty, valid from the date it is delivered by GIMA.

During the guarantee period free repair and/or replacement of any defective parts due to faulty manufacture will be given, labour, postal, transport and packaging charges etc. are not included. The guarantee therefore excludes components subject to wear and tear such as parts in rubber or PVC, etc.

No compensation will be given for loss of use of the product.

Furthermore replacement or repairs effected during the guarantee period do not lengthen the duration of the guarantee.



This guarantee is void in the event of: repairs effected by unauthorised persons or with spare parts not approved by GIMA, damage or defects caused by negligence, blows, abnormal use of the appliance or faulty installation.

The guarantee will be invalidated if the registration number has been removed, cancelled or altered.

Faulty appliances must be returned only to the retailer where the item was bought. Any item dispatched directly to us will be refused.