powered
by **BRAUN**

OXYJET / OXYJET PRO 1000-4000

Center Type 3724
Handle Type 3766

Deutsch	6
English	11
Français	15
Español	20
Português	24
Italiano	28
Nederlands	32
Dansk / Norsk	37
Svenska	42
Suomi	46
Ελληνικά	50

Internet:

www.oralb-blendamed.de
www.oralb.com
www.braun.com
www.service.braun.com

DE	AT	00 800 27 28 64 63
CH		08 44 - 88 40 10
GB		0 800 731 1792
IE		1 800 509 448
FR		0 800 944 803 (service et appel gratuits)
BE		0 800 14 592
LU		800 21172
ES		900 814 208
PT		808 20 00 33
IT		800 124 600
NL		0800-4 45 53 88
DK		70 15 00 13
NO		22 63 00 93
SE		020 - 21 33 21
FI		020 377 877
GR		800 801 3457

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English

Welcome to Oral-B!

Please read the usage instructions carefully before first use and save this manual for future reference.

IMPORTANT

- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- The toothbrush is not intended for use by children under age three. Children aged from 3 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the toothbrush, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be performed by children.
- The oral irrigator is not intended for use by children under age 8. Children aged from 8 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the product, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Oral irrigators shall be used only with liquids defined in the manufacturer's instructions.

- Use this product only for its intended use as described in this manual. Do not use attachments or chargers which are not recommended by the manufacturer.

WARNING

- If the product is dropped, the brush head/ irrigator jet should be replaced before the next use even if no damage is visible. Replace brush head every 3 months or sooner if brush head becomes worn.
- Small parts may present a choking hazard.
- Do not place the appliance in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not open and disassemble device. For battery recycling please dispose whole device according to local environmental regulations. Opening the handle will destroy the appliance and invalidate the warranty.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.

Description

A Irrigator

- 1 Jet compartment
- 2 Irrigator jet
- 3 Jet switch
- 4 Jet release button
- 5 On/off switch
- 6 Water pressure control dial
- 7 Handle switch (water stop and go)
- 8 Irrigator handle
- 9 Water container (600 ml)
- 10 Cord storage

B Toothbrush (depending on model)

- 11 Charger base
- 12 Low Charge Indicator light

- 13 Charge Indicator light
- 14 Toothbrush handle (rechargeable)
- 15 On/off button (mode selection)
- 16 Pressure sensor light
- 17 Brush head
- 18 Brush head compartment

Note: Content may vary based on the model purchased.

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Noise level: 68 dB (A)

A Using your irrigator

Before using the irrigator first time, fill the water container (9) and allow the water to run through completely by holding the irrigator handle (8) over a washbasin.

- Attach one of the irrigator jets (2) to the handle until it locks into place (a). Fill the water container with lukewarm water. If you want to add mouthwash into the water, please make sure you fill the water container with lukewarm water first, then add the mouthwash (b).
- Turn on the irrigator by pressing the on/off switch (5) (c). The handle switch should be in the lower position = water stop.
- Select the water pressure with the pressure control dial (6) (1 = gentle, 5 = strong) (d). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
- Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch (7) up to start the water stream (e). Keep your mouth slightly open to allow the water to run into the washbasin.

You can choose between two different jet functions by sliding the jet switch (3) down or up (before doing so, turn off the water stream by sliding the handle switch down):

Rotating stream for general cleaning and gum massage (jet switch in lower position): The rotation of the water stream and its micro-bubbles is designed for easy penetration and gentle gum massage (g).

Mono stream (jet switch in upper position): The straight water stream is designed for targeted cleaning (h).

Maximal operation time: 15 minutes

Cooling off time: 2 hours

In general we recommend to replace the irrigator jet every 3 months.

After use

Always empty the water container entirely. If your dentist has recommended the use of mouthwash with your irrigator, allow some clean water to run through the irrigator before you empty the container entirely to prevent any deposit build-up. Keep the irrigator switched on until the water stops running through the irrigator jet, to prevent bacteria from collecting in the standing water. Turn off the irrigator by switching off the on/off switch once the water container is empty. Longer operation of the irrigator without water may damage the motor.

Press the on/off switch to turn off the irrigator.

To take off the irrigator jet, press the jet release button (4) (f). Always wipe the handle and the irrigator jet dry; place the handle back onto its base and store the irrigator jet in the jet compartment (1).

Cleaning recommendations

Regularly clean the housing with a damp cloth. The water container (9) and the jet compartment (1) are removable and dishwasher-safe.

B Using your toothbrush (depending on model)

Charging and Operating

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the appliance into an electrical outlet and place the toothbrush handle (14) on the charge base (11) (i).
- The green charge indicator light (13) flashes while the handle is being charged. Once it is fully charged, the light turns off. A full charge typically takes 12 hours and enables at least 2 weeks of regular brushing (twice a day, 2 minutes).
Note: After deep discharge, the charge indicator might not flash immediately; it can take up to 30 minutes.
- If the battery is running low, the red charge indicator light (12) is flashing for a few seconds when turning your toothbrush on/off and the motor is reducing its speed. Once the battery is empty, the motor will stop; it will need at least a 30-minute charge for one usage.
- You can always store the handle on the plugged-in charger base; a battery overcharge is prevented by device.
Note: Store handle at room temperature for optimal battery maintenance.
Caution: Do not expose handle to temperatures higher than 50 °C. Recommended ambient temperature for charging is 5 °C to 35 °C.

Brushing techniques

Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head (17) to your teeth before switching on the appliance (j). When brushing your teeth with one of the Oral-B **oscillating-rotating brush heads** guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface (n). When using the Oral-B «**TriZone**» **brush head** place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush.

With any brush head do not forget to brush all 3 surfaces of your teeth: chewing, outer and inner surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» mode (depending on model) (optional in combination with an Oral-B «Sensitive» brush head).

Brushing modes (depending on model)

- «**Daily Clean**» – Standard mode for daily mouth cleaning
- «**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas
- «**Whitening**» – Exceptional polishing for occasional or everyday use
- «**Gum Care**» – Gentle massage of gums

How to switch between modes:

Your toothbrush automatically starts in the «Daily Clean» mode. To switch to other modes, successively press the on/off button (15). Switch off the toothbrush by holding the on/off button until the motor stops.

Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush all four quadrants of your mouth equally (k). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. If the toothbrush is turned off during brushing, the elapsed brushing time will be memorised for 30 seconds. When pausing longer than 30 seconds the timer resets.

Pressure sensor

If too much pressure is applied, the pressure sensor light (16) will light up in red, notifying you to reduce pressure (l). In addition the pulsation of the brush head will stop and the oscillation movement of the

brush head will be reduced (during «Daily Clean», and «Sensitive» mode). Periodically check the operation of the pressure sensor by pressing moderately on the brush head during use.

Brush heads

Oral-B offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle. Our **oscillating-rotating brush heads** can be used for tooth-by-tooth precision cleaning.

Oral-B CrossAction brush head

Our most advanced brush head. Angled bristles for a precise clean. Lifts and powers away plaque.

Oral-B FlossAction brush head

features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

Oral-B 3D White brush head

features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B «3D White» brush head.

Oral-B Sensi Ultrathin brush head

features a combination of regular bristles to clean teeth surfaces and ultrathin bristles for a gentle brushing experience on the gumline.

Most Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums.

We do not recommend using the Oral-B «FlossAction» or Oral-B «3D White» brush head with braces. You may use the Oral-B «Ortho» brush head, specifically designed to clean around braces and wires.


Cleaning recommendations

After brushing, rinse your brush head carefully under running water for some seconds with the handle switched on. Switch off the handle and remove the brush head. Clean both parts separately under running water (m); then wipe them dry before reassembling the toothbrush. Put the toothbrush handle back on the charger base. The charger base should be cleaned with a damp cloth only.

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

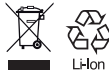Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

30 Day Risk Free Trial

Satisfaction guaranteed or your money back!
Full details at www.oralb.com or call:

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